

## Compliance Information

### FCC Compliance Statement

*This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.*

*Any modifications to this product may violate the rules of the Federal Communications Commission and make operation of the product unlawful.*

### Declaration of Conformity



Midmark Corporation declares that 9A429001 thru 9A429004 conform to the applicable requirements of the Radio Equipment & Telecommunications Equipment Directive 1999/5/EC, Annex II

Serial Number: **ALL** to present

#### Quality Assurance Manager

Don Carr  
Midmark Corporation  
60 Vista Drive  
P.O. Box 286  
Versailles, OH 45380  
U.S.A.

077-0589-00

## Limited Warranty

### SCOPE OF WARRANTY

Midmark Corporation ("Midmark") warrants to the original purchaser its new Alternate Care products and components (except for components not warranted under "Exclusions") manufactured by Midmark to be free from defects in material and workmanship under normal use and service. Midmark's obligation under this warranty is limited to the repair or replacement, at Midmark's option, of the parts or the products the defects of which are reported to Midmark within the applicable warranty period and which, upon examination by Midmark, prove to be defective.

### APPLICABLE WARRANTY PERIOD

The applicable warranty period, measured from the date of delivery to the original user, shall be one (1) year for all warranted products and components.

### EXCLUSIONS

This warranty does not cover and Midmark shall not be liable for the following: (1) repairs and replacements because of misuse, abuse, negligence, alteration, accident, freight damage, or tampering; (2) products which are not installed, used, and properly cleaned as required in the Midmark "Installation" and or "Installation / Operation Manual for this applicable product. (3) products considered to be of a consumable nature; (4) accessories or parts not manufactured by Midmark; (5) charges by anyone for adjustments, repairs, replacement parts, installation, or other work performed upon or in connection with such products which is not expressly authorized in writing in advance by Midmark.

### EXCLUSIVE REMEDY

Midmark's only obligation under this warranty is the repair or replacement of defective parts. Midmark shall not be liable for any direct, special, indirect, incidental, exemplary, or consequential damages or delay, including, but not limited to, damages for loss of profits or loss of use.

### NO AUTHORIZATION

No person or firm is authorized to create for Midmark any other obligation or liability in connection with the products.

**THIS WARRANTY IS MIDMARK'S ONLY WARRANTY AND IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED. MIDMARK MAKES NO IMPLIED WARRANTIES OF ANY KIND INCLUDING ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, THIS WARRANTY IS LIMITED TO THE REPAIR OR REPLACEMENT OF DEFECTIVE PARTS.**

## General Information

### Intended Use

*This product is intended to control patient positioning while providing better patient access by eliminating the cord from the foot control.*

### Authorized CE Representative

*Countries in the EEC should direct all questions, incidents, and complaints to Midmark's authorized CE Representative listed below:*

Promotal  
22, Rue de St. Denis  
53500 Ernee,  
FRANCE  
Phone: 33-2-430-517-76  
Fax: 33-2-430-569-10

### Calling for Service

*If service is required, contact your Midmark dealer.  
To contact Midmark directly:*

1-800-MIDMARK (1-800-643-6275)  
Monday-Friday 8:00 am until 5:00 pm (EST)  
[www.midmark.com](http://www.midmark.com)

Midmark Corporation  
60 Vista Drive  
Versailles, OH 45380-0286  
Phone: 937-526-3662  
Fax: 937-526-5542  
[www.midmark.com](http://www.midmark.com)



Because we care.



## User's Guide 625 Wireless Controls 9A429



*This guide outlines basic operation of the wireless hand / foot control and the associated base station. If service is required, contact your authorized Midmark dealer.*

**Dealer (name / phone):**

---

---

---


---

**Date of Purchase:** \_\_\_\_\_

# Association Procedure

The wireless controller will only function with the chair it has been “associated” with. This prevents unwanted movement in locations with more than one chair.

Each wireless controller’s association is preset at the factory. This procedure is only required in the event of a malfunction, if a new controller is purchased, or if you are attempting to use the controller with a chair other than the one it was originally associated with.

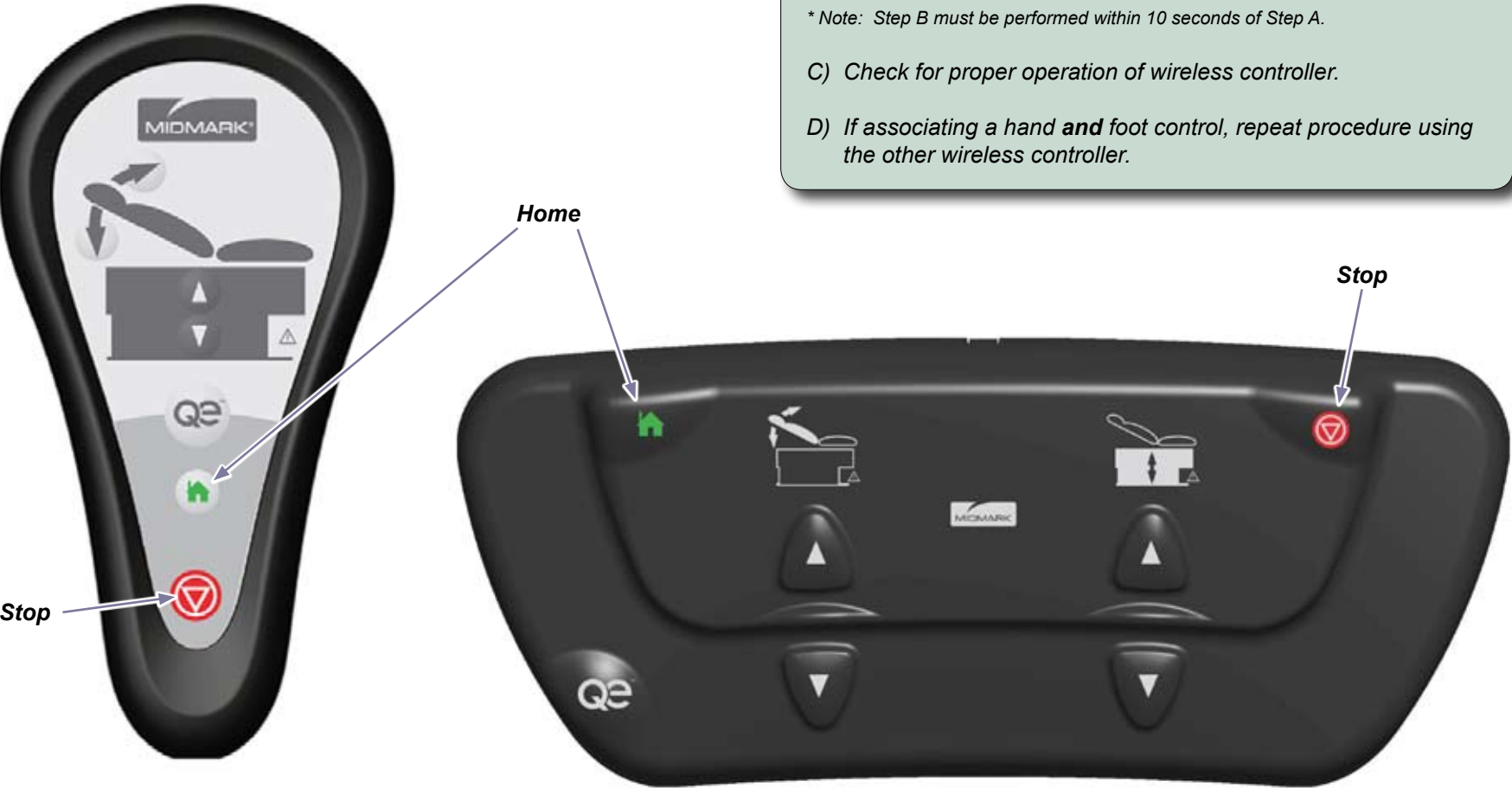


### Association Procedure:

- A) Disconnect power to the table for at least 3 seconds.
- B\*) Reconnect power, then press & hold the **Stop & Home** buttons on the **wireless** controller, until you hear three “beeps”.
- \* Note: Step B must be performed within 10 seconds of Step A.
- C) Check for proper operation of wireless controller.
- D) If associating a hand **and** foot control, repeat procedure using the other wireless controller.

# Specifications

- Range of Operation..... 20 ft. (6 m) radius
- Batteries:
- Hand Control.....AA (qty. 2)
  - Foot Control.....AA (qty. 2)

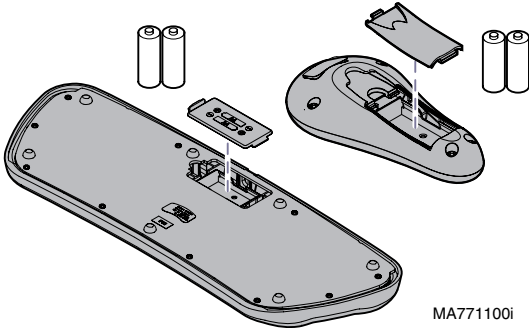

# Troubleshooting

When any control button is pressed, that controller sounds a single “beep”. (All functions operate)

- Cause: Low battery
- Solution: Replace batteries

When any control button is pressed, nothing happens. (No “beeps”, no movement, etc.)

- Cause:
- 1. Batteries completely dead
  - 2. Wireless controller has lost “association” to chair
- Solution:
- 1. Replace batteries
  - 2. Perform **Association Procedure\***

The wireless controllers will only function with the ... chair they have been “associated” with. This prevents unwanted movement in locations with more than one chair.

If further service is required, contact your authorized Midmark dealer.